COPILOT

Health Management System Version 4.1

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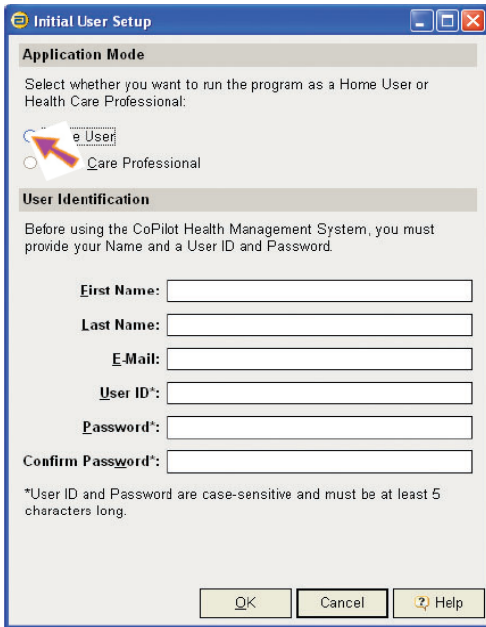
Initial User Setup

Application Mode
Select whether you want to run the program as a Home User or Health Care Professional:

Home User
 Care Professional

User Identification
Before using the CoPilot Health Management System, you must provide your Name and a User ID and Password.

First Name:
Last Name:
E-Mail:
User ID*:
Password*:
Confirm Password*:

*User ID and Password are case-sensitive and must be at least 5 characters long.

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Valitse **OK**.

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Select Unit of Measure

Select the Unit of Measure that you use on your glucose meter. Your meter is set to display test results in either milligrams per deciliter (mg/dL) or in millimoles per liter (mmol/L).

mmol/L (mmol per liter) ← Unit of Measure: mmol/L (mmol per liter) mg/dL (mg per dL)

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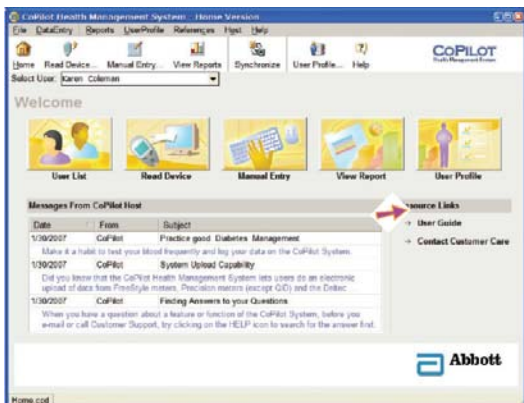
CoPilot Health Management System - Home Screen

File Data Entry Reports User Profile References Help Help

Home Read Device Manual Entry View Reports Synchronize User Profile Help

Select User: Karen Coleman

Welcome

Use List Read Device Manual Entry View Report User Profile

Messages From CoPilot Host

Date	From	Subject
1/30/2007	CoPilot	Practice good Diabetes Management
1/30/2007	CoPilot	System Update Capability
1/30/2007	CoPilot	Finding Answers to your Questions

Resource Links
→ User Guide
→ Contact Customer Care

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